

AlteCap™

Screw cap handling

Flexible systems for screw capping tasks

Automating screw cap decapping and capping is now easier and more affordable with AltemisLab's comprehensive range of decapping and capping systems. Accommodating all throughput requirements, the AlteCap™ decappers can help with the time challenges screw cap tube users often experience. All systems in the AlteCap™ range are compatible with a variety of different tube types and rack formats from any manufacturer.

Transforming lab technique.


AlteCap™ Solo



AlteCap™ Switch



AlteCap™ Swift

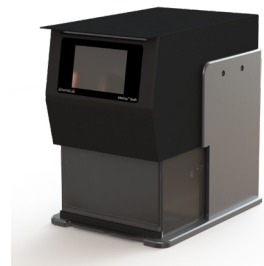

| | AlteCap™ Solo | AlteCap™ Switch | AlteCap™ Swift |
|---------------------------|---------------------|---------------------|---------------------|
| Tube throughput | 1 | 6, 8 or 12 | 24, 48 or 96 |
| Operation modes | Cap – Decap – Recap | Cap – Decap – Recap | Cap – Decap – Recap |
| Type | Handheld | Handheld | Benchtop |
| Compatible with SBS racks | Yes | 24, 48 or 96 | 24,48 or 96 |
| Tube manufacturer | Any | Any | Any |
| Fresh capping | Cap carrier | Cap carrier | Cap carrier |
| Stand included | Yes | Yes | N/A |
| Integration friendly | N/A | N/A | Yes |
| Dimensions - cm | 2.5 x 19 | 19 x 12.5 x 82 | 36.5 x 24.7 x 33.2 |
| Power | USB charging | USB charging | Mains powered |

NEW AlteCap™ Swift

Automated screw cap handling for any manufacturer's tube types

The AlteCap™ Swift decapper is a fully automated, fast & flexible capping & decapping system. Easy to change cassettes can be used with a variety of tube types, including internal & external thread tubes, from a range of different manufacturers.



Fast screw cap handling

Decap a full rack of 96 tubes in less than 20 seconds

Flexible

Easily switch between different tube & cap types

Advance cap drivers

Compatible with all screw caps from different manufacturers

Highly engineered

Individual motors apply the same torque to provide excellent seal quality

Versatile

Compatible with cap carriers for efficient capping of tubes

Compact design

Ideal for use on the bench or in the hood

Intuitive

Select from pre-set programmes to fully decap, loosen, store or discard caps, or recap tubes

Sample safety




Drip catcher prevents cross-contamination of samples during cap removal

User friendly controls

Easy to use software with touch screen display

Integration friendly

Extended rack position for easy access & remote control capability

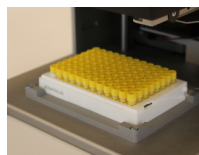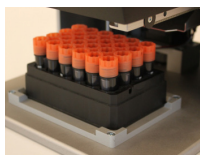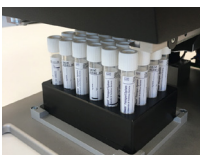

Choice of models

All of the features for 24, 48 or 96 SBS rack formats